

Product Datasheet / Instruction for use



SEC-182980 2023-02-23 1 / 2

# Rabbit IgG anti-Horse IgG (F(ab')2)-HRPO, MinX none

### **General Information**

Catalog Number: SEC-182980 Physical State: Lyophilized

Quantity: 1,5 mg Reconstitution Volume: 1.0 mL

**Reconstitution Buffer:** Restore with deionized water (or equivalent)

**Shipping Conditions:** at room temperature

**Product Expiration** Expiration date is one (1) year from date of

opening.

### **Antibody Host / Format**

Host Species: Rabbit Clonality / Isotype / Clone: polyclonal Ig

Antibody Format: IgG whole molecule Concentration: 1.5 mg/mL by UV absorbance at 280 nm

Conjugation: Horseradish Peroxidase (HRP)

Maximum Excitation / Emission: /

# **Specificity**

**Target Species:** Horse | Immunogen: Horse IgG F(ab')2 fragment

**Antibody Specificity:** IgG (F(ab')2 This product was prepared from monospecific antiserum by

immunoaffinity chromatography using Horse IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Peroxidase, anti-Rabbit Serum, Horse IgG, Horse IgG F(ab')2 and Horse Serum. No reaction was observed against Horse

Store vial at 4° C prior to restoration. For extended storage aliquot

contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at

undiluted liquid. Dilute only prior to immediate use.

room temperature. This product is stable for several weeks at 4° C as an

IgG F(c).

Minimal Cross Reactivity: none

# Formulation, Transport and Storage

**Storage Buffer:** 0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2

**Stabilizer:** 10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease

Fragment)

free

**Preservatives:** 0.01% (w/v) Gentamicin

Sulfate

Storage Temperature: 4-8°C

# **Application Recommendation**

This product is for In Vitro experimental use only. Not for Therapeutic or Diagnostic use.

Manufacturer

BIOZOL Diagnostica Vertrieb GmbH Leipziger Straße 4 85386 Eching Phone +49 (89) 3799 666 6 Fax +49 (89) 3799 666 99 E-mail info@biozol.de www.biozol.de



Your Secondary Antibody Portal Powered by BIOZOL








www.tuv.com ID 9000019771



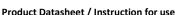



SEC-182980 2023-02-23 2 / 2

Application: ELISA, Immunohistochemistry (IHC), Immunohistochemistry (Paraffin-embedded Sections),

Immunohistochemistry (frozen sections), Western Blot

**Dilution:** ELISA 1:10.000 - 1:50.000, Immunohistochemistry (IHC) 1:500 - 1:2,500, Western Blot (WB) 1:1.000 - 1:5.000

### **Application Note:**

Anti-Horse IgG F(ab')2 Peroxidase conjugate is suitable for immunoblotting (western or dot blot), ELISA, immunoelectron microscopy and immunohistochemistry as well as other antibody-based enzymatic assays requiring lot-to-lot consistency.

# **Background Information**

Anti-Horse IgG F(ab')2 Peroxidase Antibody generated in rabbit is a proteolytic fragment of immunoglobulin G (IgG) obtained by limited digestion with the enzyme pepsin under controlled conditions of temperature, time and pH. F(ab')2 molecules lack the Fc portion of IgG and therefore receptors that bind horse IgG F(c) will not bind horse IgG F(ab')2 molecules. Secondary Antibodies are available in a variety of formats and conjugate types. When choosing a secondary antibody product, consideration must be given to species and immunoglobulin specificity, conjugate type, fragment and chain specificity, level of cross-reactivity, and host-species source and fragment composition.

### Safety Information

This reagent contains Gentamicin Sulfate (a poisonous and hazardous substance) as preservative. Although this concentration is not regarded as dangerous to health, appropriate measures should be taken to avoid skin and eye contact, inhalation, and ingestion. Standard Laboratory Practices should be followed when handling this material.DO NOT ADD SODIUM AZIDE! Use of Sodium Azide as a preservative will substantially inhibit the enzyme activity of horseradish peroxidase!

# Disclaimer

This product is warranted to perform in conformance with product specifications and to be free from defects in material and workmanship. Products are supplied for research use only. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated, this warranty is limited to one year from date of shipment when the product is subjected to normal, proper and intended usage. Buyer's exclusive remedy for non-conforming products during the warranty period is limited to replacement of or refund for the non-conforming product(s). There is no obligation to replace products as the result of accident or disaster or inappropriate use of the products in a manner for which they were not designed, or improper storage and handling of the products.







